



AUTO SEALER

REF : AS-01

PORTABLE SEALER / More than **1200** sealing cycles




Battery charge indicator



OPTION Elastomer ring avoiding slides



Allow sealing close to the blood bag

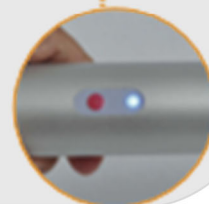

Visual LED controls the sealing cycle

Ergonomic design for easy cleaning

- Easy to use Constant sealing
- Low power consumption
- Light weight (less than 1,25 kg)
- Constant sealing quality
- level of batterie charge controled 6 LED
- Able to de used during the blood donation process




High quality of sealing and easy separation

+390544 50 161 7

export@vasinistrumenti.it


VASINI STRUMENTI

- Total control of sealing process by Micro-Processor
- Ecosystem friendly "Lithium Ion Phosphate" battery and high capacity for safety and extended lifetime

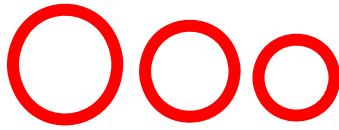

Seal all kind of PVC and EVA blood bag tubes (up to 8 mm) with motorized clamping system

TECHNICAL FEATURES

Main power: AC/DC charger Input: 100-240 VAC, 47-63 Hz, 0,3 - 0,9A Max Output: 12 Vdc / 5.25 A max Battery: Rechargeable battery 14.8 V - 3300 mA Lithium-ION Latest technology with nanophosphate materials 80 W - 50 ohm RF output max: RF frequency: 40,685 MHz Operating temperature: 10°C to 45°C (50°F à 113°F) Humidity: 10% - 90% (without condensation) Storage temperature: -10°C to 50°		
BASE (Aluminium) Dimensions : 21,7 x 8,7 x 4 cm Weight 0,75 kg	CARRYING CASE Dimensions: 44,5 x 37 x 11 cm Weight: 1,20 kg	
CLAMP (Aluminium) Dimensions : 23,3 cm (long), Ø 3,5 mm Weight 0,33 kg	PACKING Dimensions: 47 x 39 x 14 cm 2,54 kg (with case) Weight: 1,23 kg (without case)	
Sealing time: Less than 2 seconds		
CERTIFICATIONS: ISO 9001, ISO 13845, Calibration and electrical certificate CLASS I TYPE B		

ACCESORIES

		
Transport suitcase	Suitcase for 3units	Battery

Blood Bank Supplies

Via Faentina n. 207/Q – 48124 RAVENNA – ITALY T: +39
 0544 46404 F: +39 0544 501617
www.vasinistrumenti.it export@vasinistrumenti.it

